

AIRCAST®

10° Wedge Kit

LEVEL UP YOUR WALKER THE EASY WAY



enovis™

10° Wedge Kit



Level up your walker the easy way

This set of 10° heel wedges allows the Aircast® AirSelect™ pneumatic walker boot to be adapted easily for use in Achilles tendon rehabilitation. By inserting up to three wedges into the boot, the foot can be placed into 10° increments of plantar flexion, thereby reducing strain on the Achilles tendon during the various stages of recovery.


INTENDED USE/INDICATIONS

The Aircast 10° wedges are designed:

- to comfortably position the ankle joint in plantar flexion during Achilles tendon rehabilitation
- for use with the Aircast AirSelect Standard and Elite walker boots

FEATURES AND BENEFITS

- **Customized height adjustments** - the kit contains 3 wedges:
 - 3 wedges = angle of 30° of plantar flexion
 - 2 wedges = angle of 20° of plantar flexion
 - 1 wedge = angle of 10° of plantar flexion
- **Lightweight but strong** hollow design to keep additional weight to a minimum
- Double-sided adhesive tape included with each wedge for **easy application and secure fixation**
- **Available in 2 sizes:** XS-S-M, and L-XL (to fit AirSelect Standard and Elite walkers)
- **Bilateral:** fits either foot
- **Latex free**
- Sustainable design made from **recycled material**

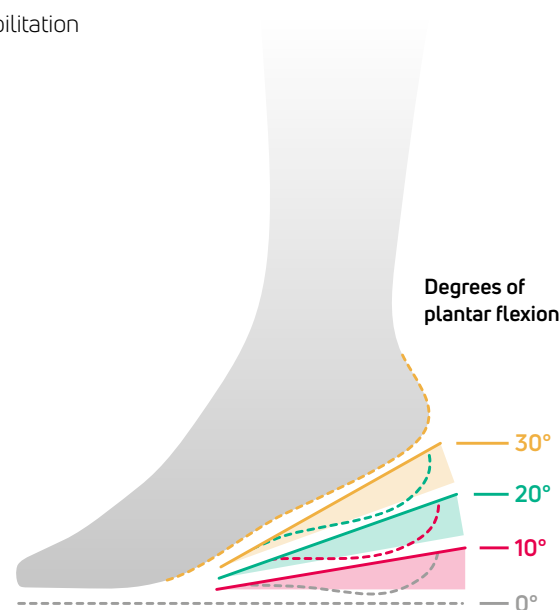

ORDERING INFORMATION

PART NUMBER	AIRSELECT WALKER BOOT SIZE
82-0110	XS-S-M
82-0210	L-XL